


# SynthoGraft™

Pure Phase Beta-Tricalcium Phosphate

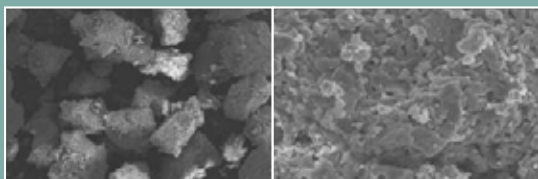


**PROCEDURE MANUAL AND CATALOG**

## Why SynthoGraft?

**SynthoGraft** offers a unique structure which provides stability, while its micro-porosity allows for rapid vascularization and subsequent resorption. Although several varieties of beta-tricalcium phosphate are now commercially available, their bone regenerating capabilities are not equal. The differences can affect not only the rate and quality of bone regeneration, but also the rate of resorption and replacement with autogenous bone during the healing process.

# SynthoGraft™

Pure Phase Beta-Tricalcium Phosphate

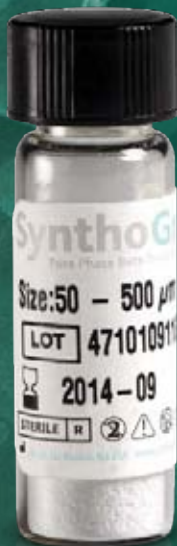

### SynthoGraft offers:

- ▶ Increased patient acceptance
- ▶ Elimination of the inherent risks associated with biologically-derived bone graft materials
- ▶ Greater surface area compared to other synthetic bone grafting materials
- ▶ Rapid vascularization and subsequent resorption when mixed with the patient's own blood
- ▶ Nanometer-scale porosity
- ▶ Available in two particle sizes: 50-500µm and 500-1000µm

# TABLE OF CONTENTS

Basic Procedures .....	4
Clinical Applications .....	5
Small Bone Defect .....	6
Sandwich Technique .....	7
Lateral Sinus Lift .....	8
Indications for Use .....	9
Warnings .....	9
Contraindications .....	9
SynthoGraft Product Information .....	10
Collagen Product Information .....	11

SynthoGraft is manufactured, packaged, sterilized and sold at ISO certified facilities.


# BASIC PROCEDURES

## PRODUCT PREPARATION:

1. Remove glass vial of SynthoGraft from sterile package.
2. Use anesthesia without epinephrine. Place patient's blood into dappen dish. Blood is most easily obtained from the surgical site or venipuncture of the antecubital fossa.
3. Progressively pour SynthoGraft into dappen dish allowing the blood to be absorbed by the material.

***SynthoGraft must only be wetted with the patient's blood.***

4. Mix SynthoGraft with a periosteal elevator for two minutes or until its consistency is putty-like.
5. SynthoGraft is ready for use when it adheres to a periosteal elevator.

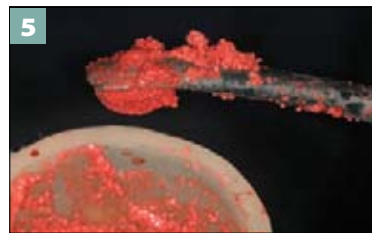


## KEYS TO SUCCESS:

- It is not recommended to mix SynthoGraft with any other bone grafting material.
- Conventional flap and curettage techniques are used to assure that the site is completely debrided and that the root surfaces are thoroughly planed and decontaminated.
- Suctioning should be limited to the excess fluid from the graft and the periphery of the site.
- Use of a membrane may be indicated.
- Clinicians are encouraged to place implants 3 to 5 months after grafting.
- Peel and use vial label for easy charting and record keeping.

## CLINICAL APPLICATIONS

### SIMULTANEOUS BUCCAL DEFECT AND INTERNAL SINUS LIFT

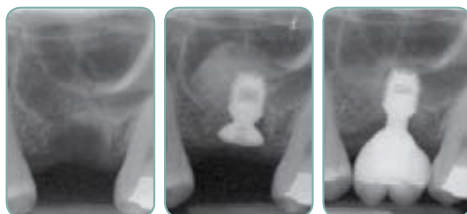

Pre-Operative

Placement

One Year

### EXTRACTION SITE




Failed Root Canal

Extraction

Post Graft

### INTERNAL SINUS LIFT

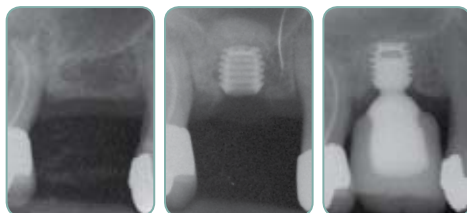


Pre-Operative

Placement

Three Years


### CREST AUGMENTATION

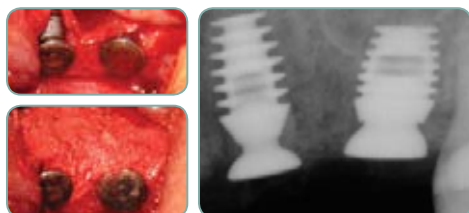

Graft In Place

Post Graft

### INTERNAL SINUS LIFT

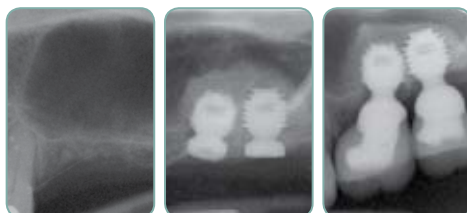

Pre-Operative

Placement

Four Years


### PERIODONTAL DEFECT



Site Of Defect

Graft In Place

### LATERAL SINUS LIFT

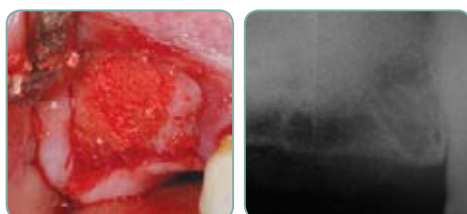


Graft In Place

Post Graft


CLINICAL APPLICATION	PARTICLE SIZE (µm)	TYPICAL QUANTITY
Extraction Site	50–500	0.25–1.00 grams
Crest Augmentation	50–500	0.25–1.50 grams
Periodontal Defect	50–500	0.10–0.50 grams
Buccal Defect/Internal Sinus Lift	50–500	0.25–0.50 grams
Internal Sinus Lift	50–500	0.25–1.50 grams
Lateral Sinus Lift	50–500 or 500–1000	1.50–5.00 grams




# SMALL BONE DEFECT

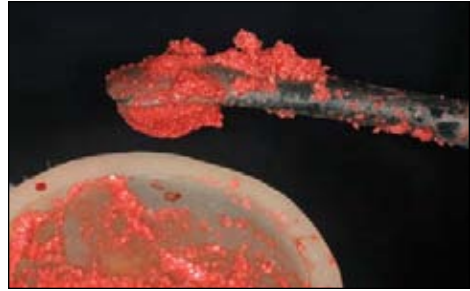

SynthoGraft is ready for use when it adheres to a periosteal elevator.

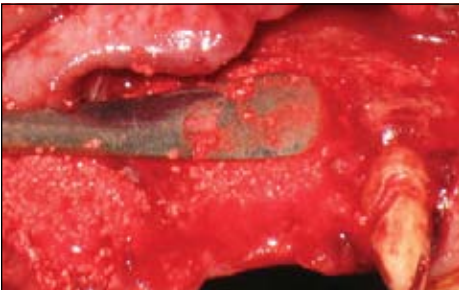

Apply SynthoGraft over the bony defect with a periosteal elevator.


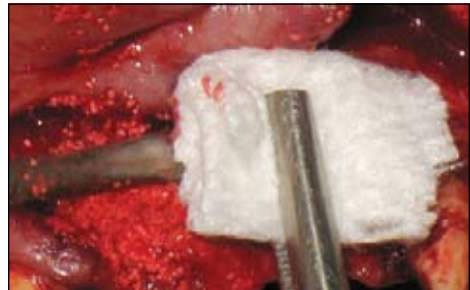


Compress SynthoGraft with periosteal elevator. Use gauze to absorb excess fluid.




Compression provides mechanical stability to the graft.

# SANDWICH TECHNIQUE

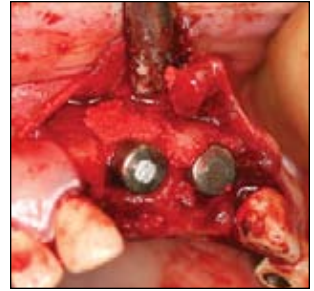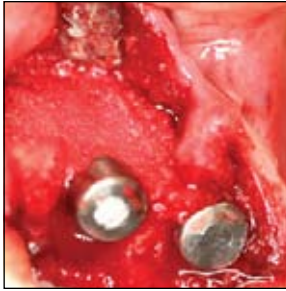

If autogenous bone is collected, place it first on the bony defect and then apply SynthoGraft on top of it. Do not mix the SynthoGraft with the harvested bone.

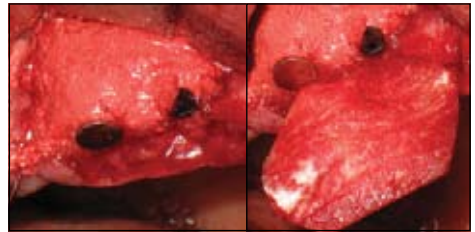

For a small-sized defect or the one-stage implant placement technique, the flap may be closed over the grafted area without a membrane.

Use a membrane for a large defect or when a horizontal augmentation of the buccal plate is required.



The mucoperiosteal flaps should be sutured to achieve primary closure.

# LATERAL SINUS LIFT



Mix SynthoGraft with a periosteal elevator for two minutes or until its consistency is putty-like.



Use a Graft Delivery Instrument to place SynthoGraft against the mesial, distal, and medial walls of maxillary sinus.



Wet resorbable collagen membrane with patient's blood prior to its placement.



Place membrane over the antrostomy.



Suture site with resorbable sutures.



## INDICATIONS FOR USE:

- Filling and/or reconstruction of a traumatic or degenerative multi-walled bone defect.
- Augmentation of the sinus floor.
- Augmentation of an atrophied alveolar ridge.
- Filling of a periodontal or other alveolar bone defect, tooth sockets, and osteotomies.
- Preservation of the alveolus for an implant osteotomy.

## WARNINGS:

- SynthoGraft should only be used by or under the supervision of trained personnel with experience in surgical techniques.
- Do not use SynthoGraft if package has been opened, damaged, or if the expiration date has passed.
- Do not compromise blood supply to the surgical site.
- Do not apply SynthoGraft unless it is wetted with the patient's blood.
- Do not wet SynthoGraft with any solution (e.g. physiological saline, NaCl or antibiotics) other than the patient's blood.
- Do not mix SynthoGraft with any other bone grafting material.
- Do not overfill surgical site.
- Do not re-sterilize SynthoGraft.
- Discard any unused SynthoGraft particles.
- Use sutures and/or membranes to prevent migration of particles.
- SynthoGraft is packaged and sterilized for single use only.

## CONTRAINDICATIONS:

SynthoGraft should not be used for patients with the juvenile form of chronic periodontitis, uncontrolled systemic diseases, infections, endocrinopathies, coagulopathies, psychological and neurological concerns, or in any other instance where the clinician believes that surgery or the use of SynthoGraft is inappropriate.

# SYNTHOGRAFT PRODUCT INFORMATION

PARTICLE SIZE (µm)	GRAMS PER VIAL	VIALS	PART NUMBER
50–500	0.25g	5	260-400-125
50–500	0.50g	5	260-400-150
50–500	1.00g	5	260-400-151
50–500	2.00g	5	260-400-152
500–1000	0.25g	5	260-400-525
500–1000	0.50g	5	260-400-500
500–1000	1.00g	5	260-400-501
500–1000	2.00g	5	260-400-502




## Dappen Dish

DESCRIPTION	PART NUMBER
Silicone Dappen Dish	260-103-030



# COLLAGEN PRODUCT INFORMATION


## Resorbable Collagen Membrane

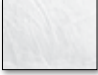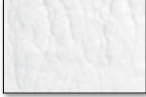DESCRIPTION	SIZE (mm)	PART NUMBER	
Resorbable Collagen Membrane	15 x 20 x 0.3	260-509-600	
Resorbable Collagen Membrane	20 x 30 x 0.3	260-509-300	
Resorbable Collagen Membrane	30 x 40 x 0.3	260-509-800	

### MEMBRANE CONSIDERATIONS


- Intended for use in oral surgical procedures as a resorbable or non-resorbable material for placement in the areas of dental implants, bone defects and ridge augmentation procedures.
- Ensure that the membrane is stable and securely placed under the tissue.
- Primary closure is desirable to allow for full coverage of the membrane.
- Dissecting and suturing the buccal aspect of the flap will allow for the coronal movement of the flap facilitating primary closure.
- Each membrane can be easily trimmed during surgery for an appropriate clinical adaptation.

## Resorbable Collagen Plug

DESCRIPTION	SIZE (mm)	PART NUMBER	
Resorbable Collagen Plug (10)	10 x 20	260-509-400	


### COLLAGEN PLUG CONSIDERATIONS

- The Resorbable Collagen Plug is fully absorbed in 10 to 14 days.
- Its cylindrical shape adapts easily to surgical sites.

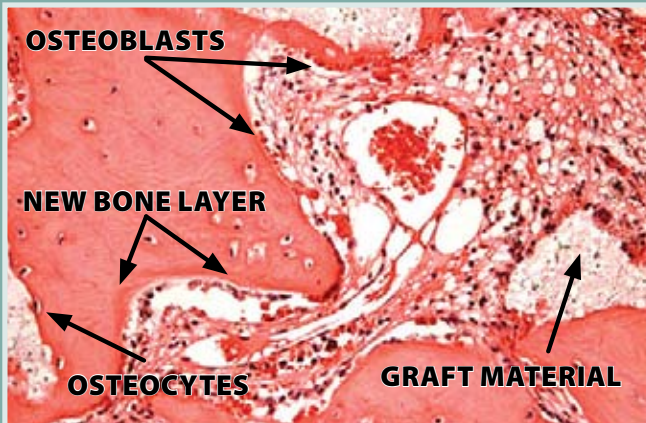

Courtesy of Ziedonis Skobe, Ph.D., Forsyth Institute and Harvard University

### 3 MONTH HISTOLOGY

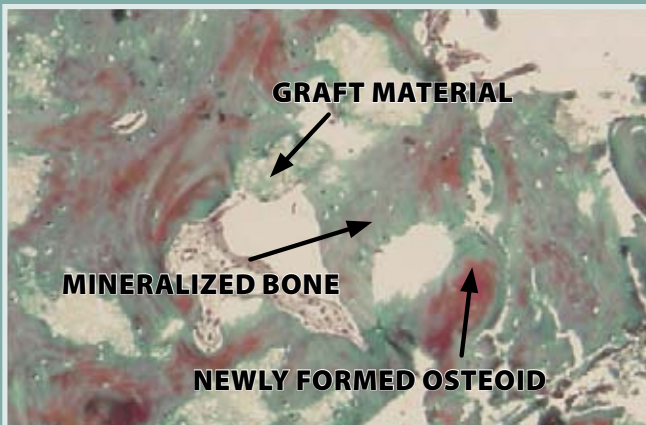

Courtesy of Serge Dibart, D.M.D., Boston University School of Dental Medicine

### 6 MONTH HISTOLOGY

**SynthoGraft**<sup>™</sup>  
 Pure Phase Beta-Tricalcium Phosphate

501 Arborway ♦ Boston, MA 02130 ♦ USA  
 TEL: 617.524.4443 ♦ FAX: 617.524.0096  
[www.synthograft.com](http://www.synthograft.com)

MADE IN THE USA